Gold Industrial Merit Award

Golden Prize of Seoul International Patent Invention Fair 2002

Golden Prize of Korea Patent Technology

Golden Prize & IFIA Special Prize of Invention & Patents in Geneva

MIGUN MEDICAL INSTRUMENT

Migun Mat Series User's Manual

MIGUN

Migun means Beautiful Health.


Please read this manual carefully for safe and effective use of this product.
Please keep this manual safe and refer to it when necessary.

Table of Contents

Thank you for purchasing a Migun Mat Series product. Please read these instructions carefully. Understanding them in whole will assist you to use this product safely and properly. It will also be helpful to keep this instruction guide handy for easy referral whenever you would need it.

*Our Migun Mat Series are listed as FDA's medical device class II. These products are intended to be used for temporary and topical muscle relaxation and relief of minor pain and stiffness.

Introduction	3
Components of the Mat	4-5
Product Specifications	6-7
Controller Parts & Function	8-11
Instructions	12-17
Warning & Caution for All Mats	18
Warning & Caution for Car-Mat	19
What are Negative Ions & Far Infrared Rays?	20-21
What is Electro-Magnetic Field?	22-23
Limited Warranty	24-25

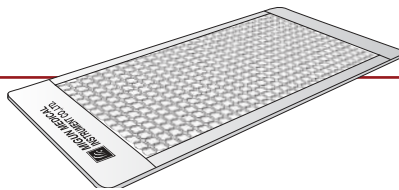MIGUN MATS

HEALTHY MAT/MINI-MAT/MINI-PAD
CAR-MAT/MINI-MAT DUO

The Migun Mats are the next steps in bringing natural health into your home. Constructed with SIX durable layers, the Migun Healthy Mat, Mini-Mat, Mini-Pad, Car-Mat and Mini-Mat Duo create healthy atmospheres for your home and your automobile by producing negative ions, generating far infrared rays, and completely neutralizing harmful electromagnetic waves. The Migun Mats are the "healthy relax system" that integrate natural and modern science that can help you attain quality undisturbed sleep and rest.

POSSIBLE BENEFITS

- Help remove cholesterol, toxins, and lipids; boost the immune system (FIR)
- Expand arteries and veins; revitalize circulation (FIR)
- Help increase metabolism and will eventually propagate cell regeneration (FIR/Negative ions)
- Successful treatment for arthritis, rheumatism and muscle spasms (FIR)
- Cells become very active; nutrition is more efficiently absorbed in cells; and the waste is effectively discharged (Negative ions)
- Allow you to achieve quality deep sleep cycles (Negative ions)

• Built-in Heating Function

• Neutralizes EMF

• Emits Negative Ions

• Emits Far Infrared Rays

*NASA concluded in the 1980's that Far Infrared stimulation is the ideal way for astronauts to maintain cardiovascular conditioning during long space flights.

**Migun Mini Pad
MG-3100**

**Migun Mini-Mat Duo
MG-3300**



**Migun Car-Mat
MG-3400**

**Migun Mini-Mat
MG-3600**


**Migun Healthy Mat
ED-301**

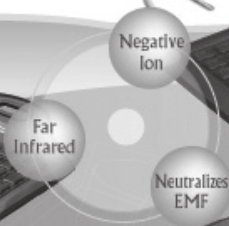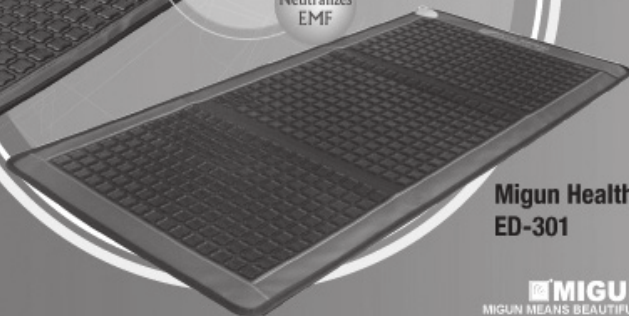
MIGUN

MIGUN MEANS BEAUTIFUL HEALTH.


Components of the Mat

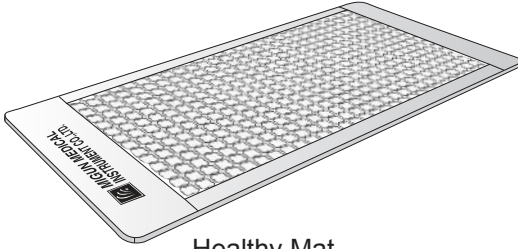Healthy Mat
(Single)

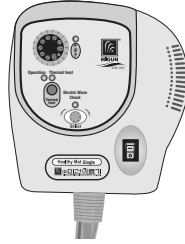Controller

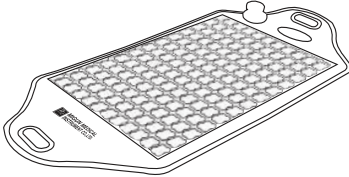Mini-Mat

Controller

Mini-Pad

Controller

Components of the Mat

Car-Mat

Controller

Mini-Mat Duo

Controller

Car Adapter (For Car-Mat only)

*Home Adapter (not included)

Manual

Product Specifications

Model Number ED-301 (Healthy Mat Single), MG-3600 (Mini-Mat)

Product Name Migun Healthy Mat, Migun Mini-Mat

Voltage/Frequency 120V, 60Hz

Power Consumption ED-301: 240W, MG-3600: 120W

Heating Emitting Wire Non-magnetic thermal wire


Temperature Adjustment Method Thermostat

Dimension/Weight ED-301: 30 lbs / 34"x75", MG-3600: 11.9 lbs / 20.5"x42.5"

Mat Temperature 86°F~140°F

Others

1. Built-in circuit to prevent excessive temperature inside the mat (using bimetal 205°F).
2. Use of natural fiber and copper mesh netting to block earthwave and electromagnetic waves.
3. Use non-magnetic thermal wire that does not emit magnetic fields

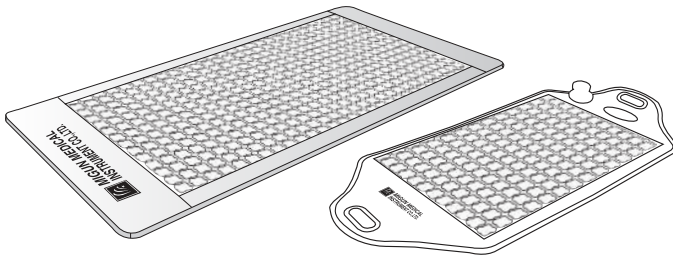
Product Specifications

Model Number	MG-3100 (Mini-Pad), MG-3400 (Car-Mat), MG-3300 (Mini-Mat Duo)
Product Name	Migun Mini-Pad, Migun Car-Mat, Migun Mini-Mat Duo
Voltage/Frequency	120V, 60Hz, 12V for Car-Mat
Power Consumption	MG-3100: 60W, MG-3400: 30W, MG-3300: 85W
Heating Emitting Wire	Non-magnetic thermal wire
Temperature Adjustment Method	Thermostat
Dimension/Weight	MG-3100: 2 lbs / 9"x25", MG-3400: 4 lbs / 17.5"x26" MG-3300: 9.9 lbs / 18.5"x18.5"
Mat Temperature	MG-3100: 86°F~140°F, MG-3400: 86°F~122°F, MG-3300: 86°F~158°F
Others	<ol style="list-style-type: none"> 1. Built-in circuit to prevent excessive temperature inside the mat (using bimetal 205°F). 2. Use of natural fiber and copper mesh netting to block earthwave and electromagnetic waves. 3. Use non-magnetic thermal wire that does not emit magnetic fields


Controller Parts & Function(Healthy Mat)

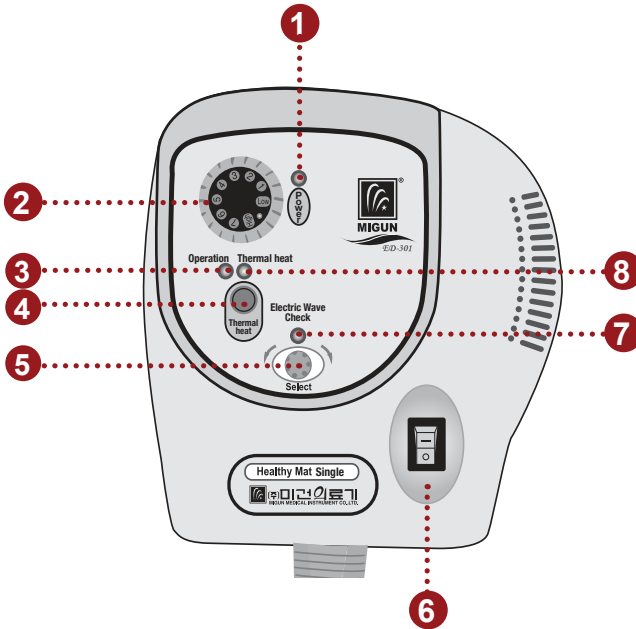
- 1. Power Indication Lamp:** Lights up when the main power is on.
- 2. Temperature Dial:** Sets temperature of the the mat.
- 3. Operation Lamp:** Lights up when the mat is heated.
- 4. Thermal Heat Start Button:** Turns on heat on the mat.
- 5. Electric Wave Select Dial:** Turn this dial to block off electromagnetic wave.
- 6. Main Power Switch:** Turns main power on or off.
- 7. EMF Indication Lamp:** This light turns off when EMF is being blocked.
- 8. Thermal Heat Lamp:** Lights up when the thermal heat start button is pressed.

Controller Parts & Function(Mini-Mat)

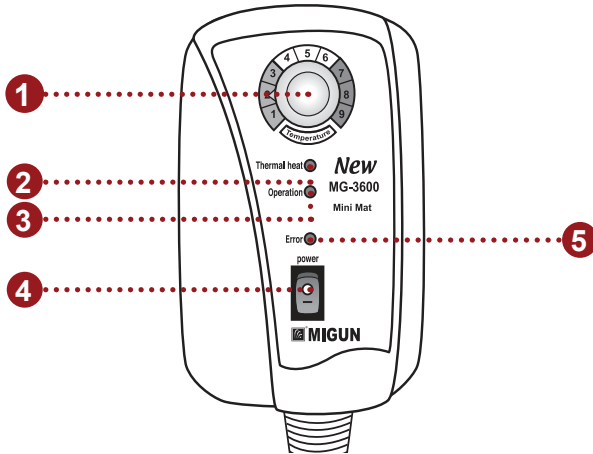1. **Temperature Control Dial:** Sets the temperature of the mat.
2. **Thermal Heat Lamp:** Lights up when power is turned on.
3. **Mat Operation Lamp:** Lights up when the mat is being heated.
4. **Main Power Switch:** Turns main power on or off.
5. **Error Lamp:** Lights up to indicate an error.

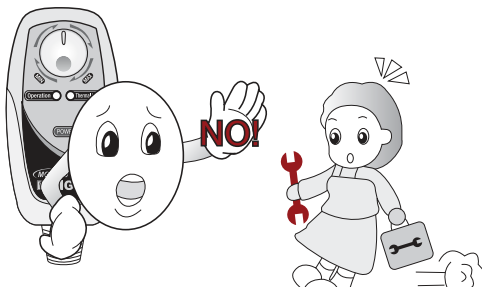*Never open or disassemble the controller. This may cause electric shock.

Controller Parts & Function(Mini-Pad/Car-Mat)

Mini-Pad

- 1. Temperature Control Dial:**
Controls temperature of the mat.
- 2. Mat Operation Lamp:**
Lights up when the mat is being heated.
- 3. Thermal Heat Lamp:**
Lights up when power is on.
- 4. Main Power Switch:**
Turns power on/off for the mat.


Car-Mat

- 1. Temperature Control Dial:**
Controls temperature of the mat.
- 2. Thermal Heat Lamp:**
Lights up when the power is on.
- 3. Mat Operation Lamp:**
Lights up when the mat is being heated.
- 4. Main Power Button:**
Turns power on/off for the mat.

Controller Parts & Function

Controller Parts & Function(Mini-Mat Duo)

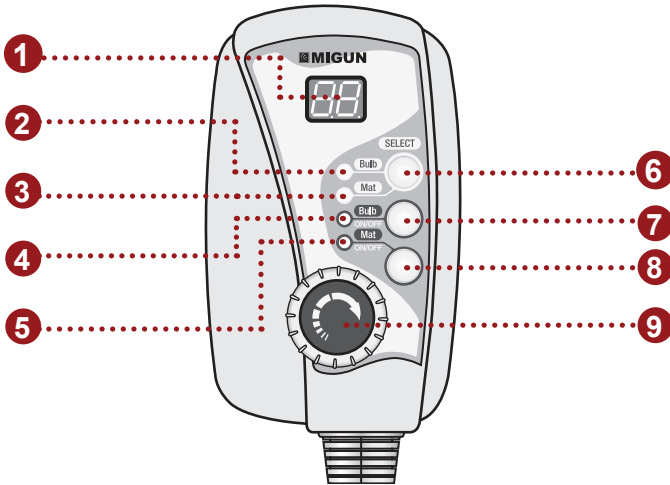
1. **Display Monitor:** Indicates temperature and operational status.
2. **Jade Bulb Lamp:** Indicates temperature of the jade bulb is being displayed on the display monitor. (Green light)
3. **Mat Lamp:** Indicates temperature of the mat is being displayed on the display monitor. (Yellow light)
4. **Jade Bulb Operation Lamp:** Lights up when jade bulb is being used.
5. **Mat Operation Lamp:** Lights up when the mat is being heated.
6. **Display Select Button:** Selects temperature for either the jade bulb or the mat to be displayed.
7. **Jade Bulb Operation Button:** Turns power on/off for the jade bulb.
8. **Mat Operation Button:** Turns power on/off for the mat.
9. **Temperature Control Dial:** Controls temperature of the mat.


Mini-Mat Duo

Instructions(Healthy Mat)

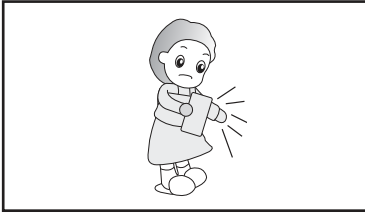
Please read the instruction guide carefully.

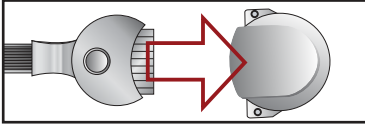
Before you connect the power plug, spread the Migun Healthy Mat on a bed or on the floor in order to smooth out any folds.

Insert the 6-pin connector of the controller into the mat.

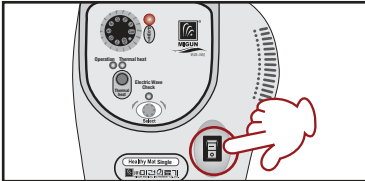
Connect the power plug into an outlet.
*Please connect to AC 120V/60Hz only.

Press the ON/OFF button. The power lamp will light up.

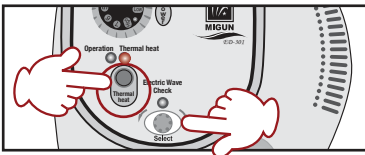
Press the thermal heat button to turn on the heat. Thermal Heat Lamp will light up. The Operation Lamp will automatically turn on or off depending on the temperature control.

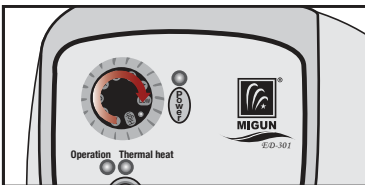
EMF indication light will turn on when harmful EMF is detected. In that case, turn the Electric Wave Select dial to turn EMF indication light OFF.

Turn Temperature Control Dial to set desired temperature on the mat. Once the product is sufficiently warm, adjust the setting for your comfort.

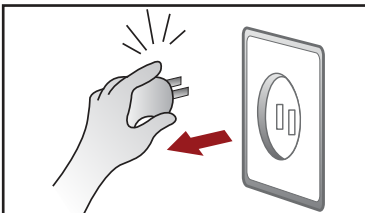
Turn Power switch off and unplug power cord from the outlet when you are not using the mat.

Instructions(Mini-Mat)

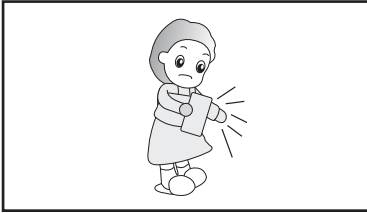
Please read the instruction guide carefully.

Before you connect the power plug, spread the Migun Mini-Mat on a bed or the floor in order to smooth out any folds.

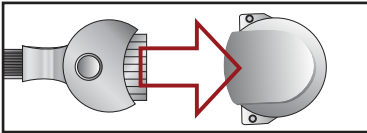
Insert the 5-pin connector of the controller into the mat.


Connect the power plug into an outlet.
*Please connect to AC 120V/60Hz only.

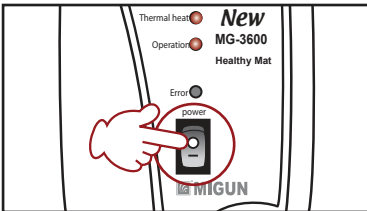Press power switch to turn on the power.
Thermal Heat Lamp will light up. The Operation Lamp will automatically turn on or off depending on the temperature control.

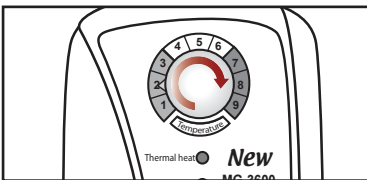
Turn Temperature Control Dial to set desired temperature of the mat. Once the product is sufficiently warm, adjust the setting for your comfort.

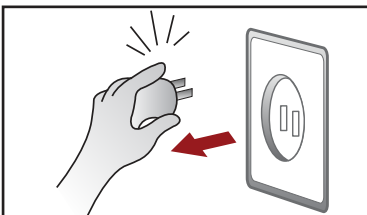
Turn Power switch off and unplug power cord from the outlet when you are not using the mat.

Instructions(Mini-Pad)

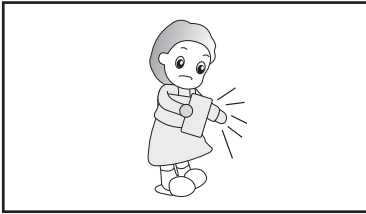
Please read the instruction guide carefully.

Before you connect the power plug, spread the Migun Mini-Pad on the body to smooth out any folds.

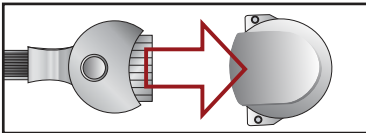
Insert the 5-pin connector of the controller into the mat.

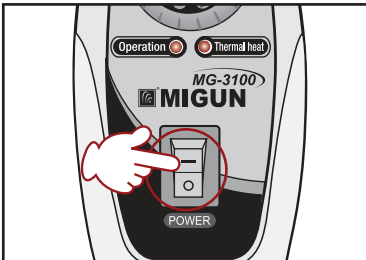
Connect the power plug into an outlet.
*Please connect to AC 120V/60Hz only.

Press power switch to turn on the power. Turn Temperature Control Dial, and Thermal Heat Lamp will light up. The Operation Lamp will automatically turn on or off depending on the temperature control.

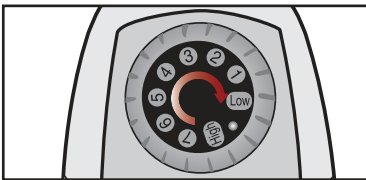
Turn Temperature Control Dial to set desired temperature of the mat. Once the product is sufficiently warm, adjust the setting for your comfort.

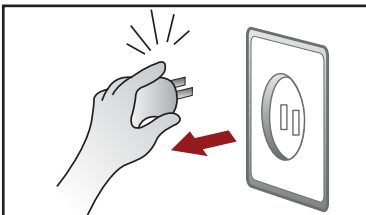
Turn Power switch off and unplug power cord from the outlet when you are not using the mat.

Instructions(Car-Mat)

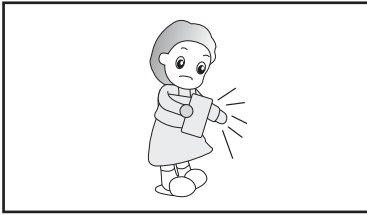Please read the instruction guide carefully.


Before you connect the power plug, spread the Migun Car-Mat on a flat surface in order to smooth out any folds.

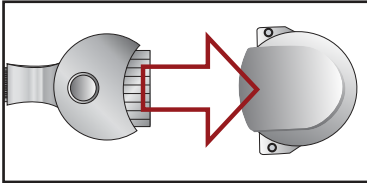Insert the connector to the adapter.

Connect the adapter into an outlet.


*This product uses 12V DC for automobiles and AC120V/60Hz for indoors.

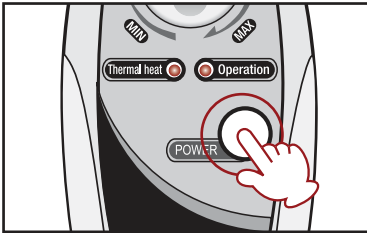Press power button to turn on the power. Turn Temperature Control Dial, and Thermal Heat Lamp will light up. The Operation Lamp will automatically turn on or off depending on the temperature control.


Turn Temperature Control Dial to set desired temperature of the mat. Once the product is sufficiently warm, adjust the setting for your comfort.

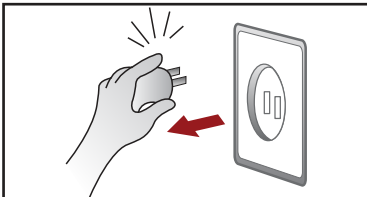
Turn Power switch off and unplug the adapter from the outlet when you are not using the mat.

Instructions(Mini-Mat Duo)


Preparation Before Use

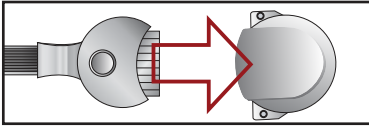
Please read the instruction guide carefully.
Insert the 7 pin connector of the controller into the mat.

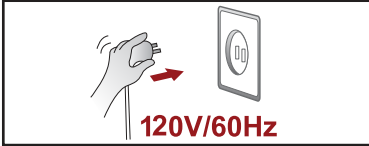Connect the power plug into an outlet.
*Please connect to AC 120V/60Hz only.

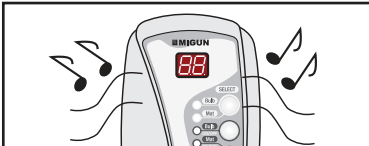
The controller will make a chime bell sound, and "oF" will show with green light on the display monitor.

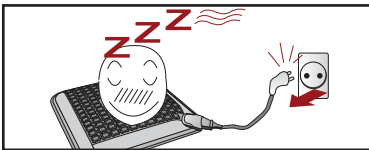*When the mat is turned on for more than 12 hours, the power will turn off automatically. Please plug off the mat and plug it back in to use MiniMat-Duo again.


How to To Use Jade Bulbs

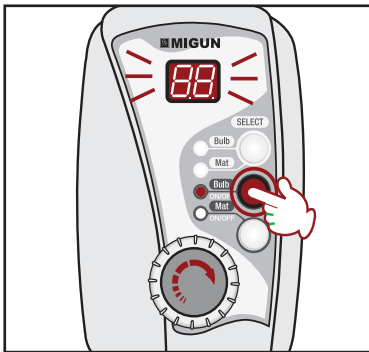1. Press Jade Bulb Operation Button to start.
 - Green Bulb Lamp will turn on, and the temperature on the display monitor will blink three times before starting.
2. Turn the Temperature Control Dial to set the desired temperature.
 - Temperature will rise between 30°C~70°C (86°F~158°F). Temperature below 30°C (86°F) will not be indicated on the display monitor.


30°C: 86°F / 40°C: 104°F / 50°C: 122°F /
60°C: 140°F / 70°C: 158°F

Instructions(Mini-Mat Duo)

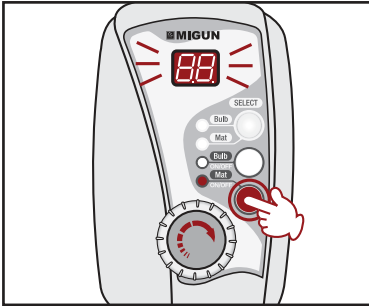How to Use the Mat

1. Press Mat Operation Button to start.
 - Yellow Bulb Lamp will turn on, and the temperature on the display monitor will blink three times before starting.
2. Turn the Temperature Control Dial to set the desired temperature.
 - Temperature will rise between 30°C~70°C (86°F~158°F).

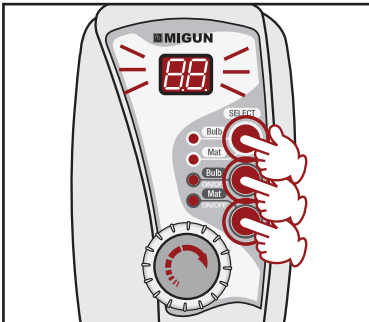How to Use Both Mat & Jade Bulbs


1. When Jade Bulb and Mat are used simultaneously, they will be controlled at the same temperature. You can set the temperature for Mat or Jade Bulbs, and it will control temperature for both.
 - You can change the temperature setting during the operation.
 - You can turn on/off the Mat /Jade Bulb during the operation.

Storage & Maintenance for All Mats

- Keep the product away from direct sunlight. Do not allow dust and foreign substances to come in contact with the product.
- Do not soak the product in water. Store the product in a dry area when not in use.
- Do not adjust or alter any components of the product under any condition.
- Please unplug the power plug when not in use and also prior to cleaning or repairing the product.
- Keeping the product clean will add longevity and continued performance to your heating mat.

Warning & Caution for All Mats

Please follow these instructions in order to protect yourself from accidents.

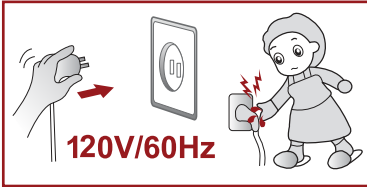
Any power input other than 120V/60Hz may result in malfunction or damage.

Please monitor the heat carefully while using the product for safety.


Please do not use a thick padding while using the product.

Do not have sharp objects puncture or damage the product.

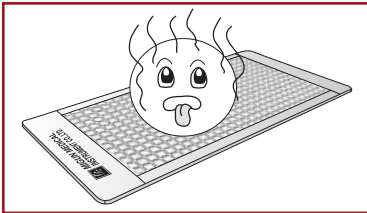Do not pull the electric cord with wet hand. This may cause electric shock.


Please make sure the product is dry before using. If you spill liquid on the product, please dry the product completely before using.

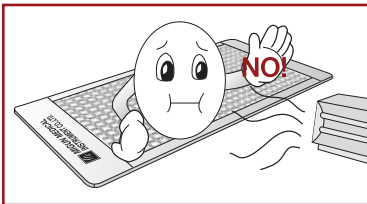Heat can cause a skin irritation or burn. Please lower the temperature if it becomes too hot.

Do not use the product folded or crumpled.

Do not use the product near any fire, heating appliances, or flammable objects. This may set fire or electric shock.

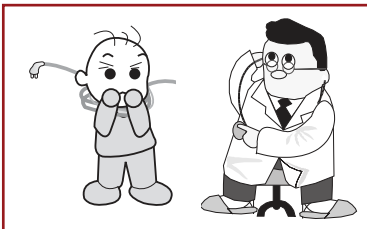Please do not let children play on or near the product.

Please stop using the product and contact Migun representative if there is any problem. Please consult with your doctor if you have any medical question.

Warning & Caution for Car-Mat

WARNING

1. Use this unit only for its intended use as described in this manual.
2. This unit should never be left unattended when plugged in. Unplug from the power source when not in use and before attaching or removing adopter.
3. Close supervision is necessary when this unit is used by or near children or disabled person.
4. Never use this unit if it has damaged cord or plug; if it is not working properly; or if it has been exposed to sunlight or if it is wet.
5. Keep the product away from heated surface or water.
6. Never drop or insert any object into any opening.
7. Do not use where aerosol spray products are being used or where oxygen level is high.
8. DO NOT operate under a blanket or pillow. Excessive heating can occur causing fire, electric shock or injury to persons.
9. Do not carry this unit by the cord or use the cord as a handle.
10. To disconnect, turn all controllers off, then remove the plug from the power source.
11. Any individual who is pregnant, has pacemaker, suffered from diabetes, phlebitis, and/or thrombosis, or is at increased risk of developing blood clots due to recent surgery should consult his/her physician before using this unit.


CAUTION

1. To reduce the risk of fire, electrical shock & damage, Do not expose this unit to Rain(water), Moisture and Direct Sun Light.
3. Never cover the unit when it is in operation.
4. Extensive use at any time could lead to excessive heating that shortens the life of this unit. Should this occur, discontinue use and allow the unit to cool before operating.
5. Never use this unit directly on swollen/inflamed areas of the skin or on area of skin eruption.
6. This is a non-professional unit designed to provide soothing low heat to your contacted areas. Do not use as substitute for medical attention.
7. This unit should never be used by any individual suffering from any physical ailment that would limit the user's capacity to operate the controls.
8. Do not leave the product on for too long. This unit could make you drowsy.
9. Do not leave the unit connected to power source when user is not in the car.
10. Do not drop the product. It might damage the unit.

What is Negative Ion?

Negative Ion Comparison Chart

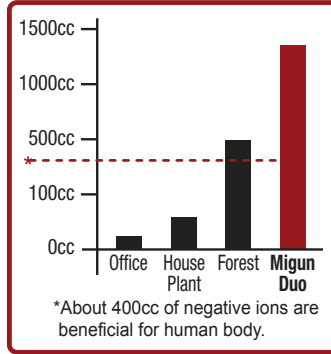
Many people think of the anion as a “mysterious” substance. However, the anion is not a substance but an electrical term that indicates a negative charge. If the molecule or atom of a substance gains one or more electron through a certain stimulation, it becomes anion or a negative ion. If it loses an electron, it becomes a cation. Therefore, there are many different kinds of negative ions according to the scientific definition.

Many researches show that oxygen negative ions enter the body through breathing, react with hemoglobin to get mixed into to blood, and are then carried to various parts of the body to facilitate metabolism. It promotes cellular activity, and as a result it boosts the immune system, among other things. When exposed frequently to an abundant supply of anion, all the bodily organs—including the nervous system, cells, blood and lymph nodes—become revitalized and strengthened. The negative ion raises the ionization ratio of the minerals in the blood, such as calcium, sodium and potassium. It also cleanses the blood through the electrical transport of substances across the cell membrane and increase the amount of globulin (the immune substance contained in the serum), thus bolstering the immune system.

In an experiment conducted in Israel, colon bacillus, protozoa and other bacteria were exposed to the atmospheric negative ions. Within six hours the number had been reduced by more than 50%. The results of a U.S. study show that when people are exposed to anionized air they lose excess weight and sleep better.

What are Far Infrared Rays?

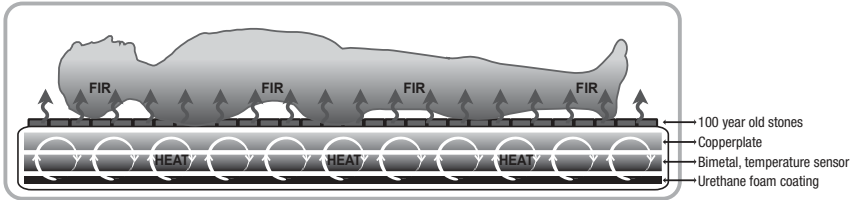Far infrared is part of the spectrum of light emanating from the sun. Far infrared is not visible to the eye, but you feel it as heat penetrating into your skin. While all light transfers energy from one place to another, infrared, literally meaning “below red”, specifies the emanation of heat energy.

When the heat system in the Migun Mat is turned on, infrared rays enter 4-5 cm into the skin. Far Infrared Rays have been proven to increase blood circulation, metabolism, and boost the immune system. They can positively affect muscles on the skin, blood vessels, lymphs, and even nervous system.


In simple terms, the Far Infrared rays cause cells in these body parts to resonant. The Far Infrared Rays heat body from within expanding arteries and veins which revitalize blood circulation, and strengthen metabolism, thus effecting chemical reactions in the cells. The results are purified blood and strengthened ability to revitalize.

*NASA concluded in the 1980's that Far Infrared stimulations is the ideal way for astronauts to maintain cardiovascular conditioning during long space flights.

What is Electro-Magnetic Field(EMF)?

Built in non-magnetic field heat rays and bimetal prevents the production of EMF.

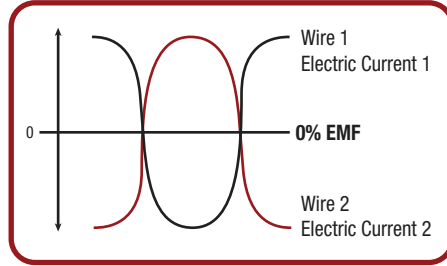Electromagnetic waves are made up of a combination of electric waves and magnetic waves. Together this is known as Electro-Magnetic Field (EMF). This EMF causes Electro-Magnetic Interference (EMI). Between the magnetic waves and electric waves, the magnetic waves have been proven to be especially harmful to your body. These are also the most difficult to neutralize in protecting our health. Studies show long-term exposure to electromagnetic waves can lead to leukemia, a lower sperm count in men, and irregular menstruation in women or possibly even birth defects. In severe cases, electromagnetic waves can cause miscarriage and brain tumors.

Migun Mat is designed to totally neutralize electromagnetic waves! When the Migun Mat neutralizes harmful EMF waves, it allows you to achieve both, the deep healing NREM (non-rapid eye movement sleep) and REM (rapid Eye Movement sleep). Both sleep stages are very important for your healing. Awaking undisturbed from REM allows you to feel more refreshed and vibrant than having your sleep cycle disrupted. Studies have also shown you age faster when you don't sleep well. So get your healthy sleep now with a Migun Mat!

Emits Negative Ions

Emits Far Infrared Rays

Protects from Electromagnetic Waves

Migun Limited Warranty (USA and Canada Only)

Written Quality Service of MIGUN Heating Mat

Product Name	Heating Mat			Model Name	
Grace Period	1 year from the date of purchase			Product Number	
Purchase Date	Month	Date	Year	Place of Purchase	
Customer Information	Name			Migun Center Name	
	Address			Address	

Migun Life, Inc. (“Migun”) warrants this product against defects in material or workmanship for the time periods and as set forth below. Pursuant to this Limited Warranty, Migun will, at its option, (i) repair the product using new or refurbished parts or (ii) replace the product with a new or refurbished product. For the purpose of this Limited Warranty, “refurbished” means a product or part that has been returned to its original specifications. In the event of a defect, these are your exclusive remedies.

Term: For a period of twelve (12) months from the original date of purchase of the product (or a period of ninety (90) days for commercial use), Migun will, at its option, repair or replace with new or refurbished product or parts, any product or parts determined to be defective.

This Limited Warranty covers only the hardware components packaged with this Product. It does not cover technical, medical, or therapeutic assistance for hardware usage.

This Limited Warranty supersedes all other warranty provided either wise including by written form or oral representations by Migun or any of its representations in relation to the purchase of the product.

Instructions: In-home and on-line diagnostic warranty service is provided during the twelve (12) month limited warranty term (ninety (90) days for commercial use) for Thermal Massage products through Migun-authorized service facility. Otherwise, to obtain warranty service, you must deliver the product, freight prepaid, in either its original packaging or packaging affording an equal degree of protection to the Migun authorized service facility specified. A dated purchase receipt or bill of sale is required. For specific instructions on how to obtain warranty service for your product, please contact any Migun-authorized service facility, or call the Migun Customer Service Center at 1-888-660-5948.

Migun Limited Warranty (USA and Canada Only)


Repair/ Replacement Warranty: This Limited Warranty shall apply to any repair, replacement part or replacement product for the remainder of the original Limited Warranty period or for ninety (90) days, whichever is longer. Any parts or product replaced under this Limited Warranty will become the property of Migun.

This Limited Warranty only covers product issues caused by defects in material or workmanship during ordinary consumer use; it does not cover product issues by any other reason, including but not limited to, product issues due to (i) acts of God, (ii) Misuse, (iii) customer instruction, installation, set up adjustments unless it was guided by Migun's on-line support, (iv) modification of or to any part of the Migun product, (v) shipping damage. This Limited Warranty does not cover Migun products sold AS IS or WITH ALL FAULTS or consumables (such as wear and tear covers). This Limited Warranty does not cover products purchased from stores other than Migun-authorized dealers. Hence, this Limited Warranty only applies when original purchase receipt signed by an authorized Migun dealer is presented. This Limited Warranty is invalid if the factory-applied serial number has been altered or removed from the product. This Limited Warranty is valid only in the United States and Canada.


LIMITATION ON DAMAGES: MIGUN SHALL NOT BE LIABLE FOR ANY INCIDENTAL OR CONSEQUENTIAL DAMAGES FOR BREACH OF ANY EXPRESS OR IMPLIED WARRANTY ON THIS PRODUCT.

DURATION OF IMPLIED WARRANTIES: EXCEPT TO THE EXTENT PROHIBITED BY APPLICABLE LAW, ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE ON THIS PRODUCT IS LIMITED IN DURATION TO THE DURATION OF THIS WARRANTY.

Some states do not allow the exclusion or limitation of incidental or consequential damages, or allow limitations on how long an implied warranty lasts, so the above limitations or exclusions may not apply to you. This Limited Warranty gives you specific legal rights and you may have other rights which vary from state to state.

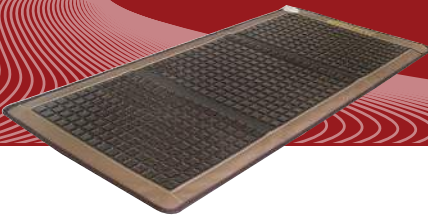Live Healthy.

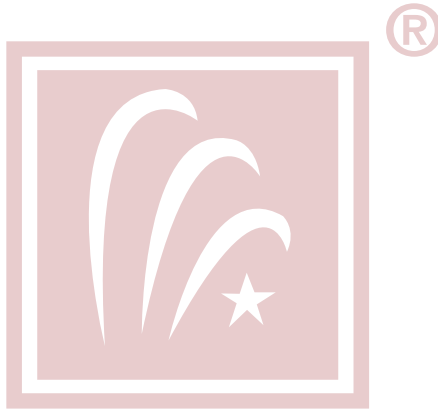 **MIGUN**
www.migunlife.com
1-828-684-1800

Gold Industrial
Merit Award

Golden Prize of Seoul
International Patent
Invention Fair 2002

Golden Prize of Korea
Patent Technology

Golden Prize & IFIA
Special Prize of
Invention & Patents
in Geneva

MIGUN MEDICAL INSTRUMENT

Migun Mat Series

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*Piedra
Cetendaria*

*No Emite
Peligroso
EMF*

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*Rayos Infrarrojos
Lejanos*

MIGUN

Migun means Beautiful Health.


Please read this manual carefully for safe and effective use of this product.
Please keep this manual safe and refer to it when necessary.

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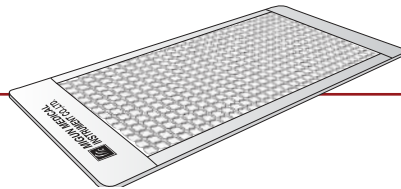
**Migun Mini-Mat
MG-3600**


**Migun Healthy Mat
ED-301**

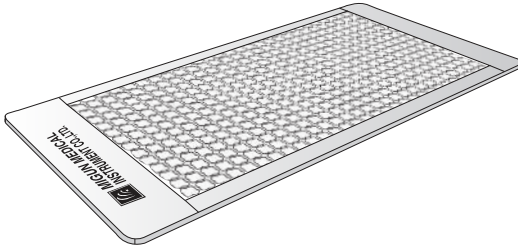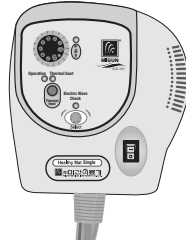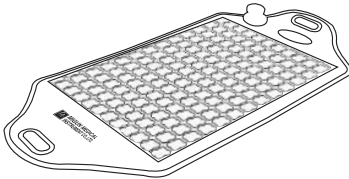
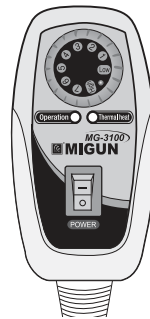
MIGUN
MIGUN MEANS BEAUTIFUL HEALTH.

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